




Sonic Travel Toothbrush

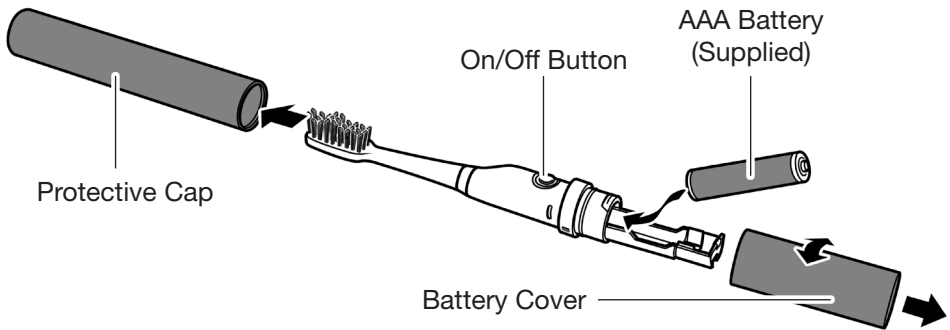

Installing The Battery

Twist the battery cover clockwise enough to release the locking thread and pull down to reveal the battery tray. Insert a 1.5V 'AAA' sized battery (supplied) matching the polar ends. Replace the battery cover and twist anti-clockwise to lock into place.

Using The Toothbrush

Remove the protective cap then apply a modest amount of toothpaste to the bristles. Press the  button to turn on the toothbrush. Move the brush head slowly and gently across the teeth in a small back and forth motion so the longer bristles can reach between your teeth. Do not scrub like a manual toothbrush. Once finished press the  button to turn the toothbrush off and rinse the bristles. Replace the protective cap.

After Sales Support

UK/N.IRELAND HELPLINE N^o: 01270 508538
REP.IRELAND HELPLINE N^o: 1800 995 036
WEBSITE: www.qesh.co.uk
EMAIL: enquiries@qesh.co.uk

Model Number: 77535

INCLUDES:
• Sonic Travel Toothbrush
• 1 x AAA battery




General Information & Safety Instructions

Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision. Children shall not play with the appliance.

This Sonic Travel Toothbrush is powered by a 'AAA' 1.5V battery.


Do not dispose of batteries in your household waste.

Please contact your local authority for details on safe disposal.

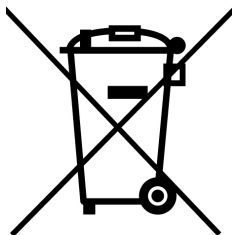

- Never throw batteries into a fire or expose to excessive heat sources.
- If batteries are swallowed, please seek immediate medical attention.
- Always ensure correct +/- battery polarity when installing batteries.
- Do not recharge non-rechargeable batteries.

Equipment, which is marked with the WEEE logo (as shown on the left), should not be thrown away with your household waste. Contact your local authority waste disposal department, as they will be able to provide details of the recycling options available in your area.

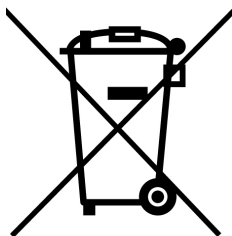

Electrical and electronic equipment (WEEE) contains materials, parts and substances, which can be dangerous to the environment and harmful to human health if the waste of electrical and electronic equipment (WEEE) is not disposed of correctly.

Equipment, which is marked with the WEEE logo (as shown on the left), should not be thrown away with your household waste. Contact your local authority waste disposal department, as they will be able to provide details of the recycling options available in your area.

After Sales Support

UK/N.IRELAND HELPLINE N°: 01270 508538
REP.IRELAND HELPLINE N°: 1800 995 036
WEBSITE: www.qesh.co.uk
EMAIL: enquiries@qesh.co.uk

Service Helpline

QESH LTD,
B7, FIRST BUSINESS PARK
FIRST AVENUE, CREWE, CHESHIRE.
CW1 6BG

VISAGE[®]

Warranty Certificate

Warranty

This product is warranted to be free from defects in workmanship and parts for a period of 36 months from the date of purchase. Defects that occur within this warranty period, under normal use and care, will be repaired, replaced or refunded at our discretion, solely at our option with no charge for parts and labour. The benefits covered by this warranty are in addition to all rights and remedies in respect of the product the consumer has under existing U.K. laws.

Proof of Purchase

This warranty is valid for the original purchase and is not transferable.

Please keep your purchase receipt as proof of purchase and as proof of the date of when the purchase was made. The receipt must be presented with the warranty card when making a claim under this warranty.

Service during the Warranty Period

Please ensure the product is properly packaged so as to ensure that no damage occurs in transit. Also make sure that you have included a detailed explanation of the problem.

Extent of the Warranty

This warranty is limited to defects in workmanship or parts. All defective products or parts will be repaired or replaced. This warranty does not extend to accessories.

Normal Wear and Tear

This warranty does not cover normal wear and tear of the product or parts.

Exclusions

This warranty does not cover the following:

- The brush head
- None oral use and use on none humans.
- Any defect caused by accident, misuse, abuse, improper installation or operation, lack of reasonable care, unauthorised modification, loss of parts, tampering or attempted repair by a person not authorised by the distributor.
- Any product that has been damaged by lightning strike either directly or indirectly, or on a main's or vehicle outlet power surge or liquid ingress.
- Any product that has not been installed, operated or maintained in accordance with the manufacturer's operating instructions provided with the product.
- Any product that has been used for purposes other than domestic use.
- The product if it is located outside of the U.K. and Ireland
- Any damage caused by improper power input or improper cable connection.

This warranty does not affect your statutory rights.

After Sales Support

UK/N.IRELAND HELPLINE N^o: 01270 508538
REP.IRELAND HELPLINE N^o: 1800 995 036
WEBSITE: www.qesh.co.uk
EMAIL: enquiries@qesh.co.uk



ViSAGE®

Sonic Travel TOOTHBRUSH

Your details:

Name _____

Address _____



_____ E mail _____

Date of purchase* _____

* We recommend you keep the receipt with this warranty card

Location of purchase _____

Description of malfunction:

Return your completed warranty card to:

Qesh Limited
B7, First Business Park
First Avenue
Crewe, Cheshire
CW1 6BG
United Kingdom
enquiries@qesh.co.uk

AFTER SALES SUPPORT

  01270 508538  www.qesh.co.uk
 1800 995 036

MODEL: 77535 PRODUCT CODE: 77535 06/2017

Phone lines available Monday to Friday, 8am - 6pm.
Calls cost 10p per minute from a landline, calls from
mobiles may vary.

3

YEAR WARRANTY